ENG

SPA

FRA

ITA

DEU

NL

POR

Shaper SH-2X

REF.: SH002X

! ENG- CAUTION

Read all precautions and instructions in this manual before using this equipment. Keep this manual for future reference.

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feel better!

www.fytter.com

EUROPEAN UNION
EUROPEAN REGIONAL
DEVELOPMENT FUND

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(ENGLISH)

Please keep and read it (safety attention) before using the caring items heredy is recored necessary important contents for safety. This notice divides into three classes when you make a mistake with handling and using.

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(ENGLISH)

Prohibited for all persons with one of these medical treatments:

- Artificial cardiac pacemaker.
- Artificial heart and/or maintained by electronic systems.
- Electrocardiograph for medical treatments.

ANY MISTAKE CAN PRESENT SIGNIFICANT DANGER. PLEASE, ASK FOR MEDICAL ADVICE BEFORE USING THE DEVICE.

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(ENGLISH)

Our bodies produce low-frequency electricity we can't feel. But electrical impulses can cause changes to the body. Let your body absorb this low-frequency electricity on targeted parts of your body to stimulate your muscles, relieve fatigue and beat stress.

- Relieve back, shoulders and waist problems.
- Improvement of muscular paralysis.
- Massage effect.
- Increase the benefits of a better diet.

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CONTENTS - CONTENIDOS - SOMMAIRE - INDICE - INHALT - INHOUD - CONTEÚDO.

(ENGLISH)

- 1.What is low frequency treatment?
2. Product
 - 2.1 Product standards
 - 2.2 Product feature
3. How to use
 - 3.1 How to operate
 - 3.2 How to use the pad
 - 3.3 Process for operation
 - 3.4 How to use and effective
 - 3.5 Taking care after using
4. How to manegment
 - 4.1 Main machine
 - 4.2 Adhesive gel pad
 - 4.3 Change the battery
 - 4.4 An effective using of adhesive gel pad

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Mode 6
○ ○ ● ●
-F-

(ENGLISH)

When you push the button "Turn on" and "Increase the intensity" (= key "ON"), the LEDs start to light up. When you push the button "Turn off" and "Decrease the intensity" (=key "OFF"), the LEDs switch off.

- Buttons ON and OFF: intensity is divided into 10 levels.
- Button PROGRAM: 6 programs (A-F).

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(ENGLISH)

How to use the pad

- 1) Clean the parts of your body that will be in contact with the device.
- 2) Pull out the gel pads from the vinyl plastic.
- 3) Put the adhesive gel pads on the targeted parts of your body (arms, shoulders ...).

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1. Don't attach it with "power on", because for shocking.
2. Don't touch attach side of adhesive pad, it make deereasing using time.

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(ENGLISH)

- 1) Take out the device from its box.
- 2) Put the adhesive gel pads on the targeted parts of your body (arms, shoulders, ...).
- 3) Turn on the device.
- 4) Choose the mode.
- 5) Adjust output level.
- 6) Working time = 12 minute sessions.
- 7) Push button OFF until the LEDs switch off.

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Please start it with at lowest output
when it is beginning

(ENGLISH) 3.4 How to use and effective.

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(FRENCH) 3.4 Comment utiliser et efficace.

(ENGLISH)

- 1) The low-frequency electricity is suitable for 12 minutes, 2 or 3 times a day.
- 2) It is not recommended to use the device in the next 30 minutes following a meal, and it is advised to use after a warm shower.

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(ENGLISH)

- 1) Switch off after use.
- 2) Separate adhesive gel pads from the device.
- 3) To avoid the gel pad drying out, stick it to the support.
- 4) The casing should be cleaned with a wet rag and a neutral detergent.
- 5) Keep the device away from sun rays, hot temperatures, humidity, dust, fire, ect ...
- 6) Keep away from children.
- 7) Don't touch the battery with wet hands.
- 8) If not used during more than a month, please remove the battery and store it separately.

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(FRENCH) 4.1 machine principale

(ENGLISH)

- 1) Please clean the device with a wet rag and neutral detergent, then dry.
- 2) Keep the device away from sun rays, hot temperature, humidity, dust, fire, ect ...

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(FRENCH) 4.2 Coussin de gel adhésif

(ENGLISH)

1. If the self-adhesive part is dirty or sticks less and less; then use a wet rag to bring the adhesiveness back.
2. Gel pads (electrodes) should be replaced by some of the same model.

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(FRENCH) 4.3 Charger la batterie

(ENGLISH)

Please connect the battery to a computer or any other USB port thanks to the USB cable.

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(SPANISH) 2.1 Normas de producto

(FRENCH) 2.1 Normes de produits

Product Name	Shaper SH-2X
Using mode	A - F (6mode)
Frequency	1~100 Hz
Working time	12 Minutes (Per installation time)
Output	Maximum 9.8 mA
Size	21 x 10 x 2,3 cm
Weight	57.4g

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(ENGLISH)

2.2 The low-frequency electricity works through rhythm and soft impulses.

The low-frequency electricity has a special rhythm and regular impulses. Depending on the power and speed, choose the pulse that adapts the best to your body. The low frequency rhythm produces a massage effect with its pressures on your body.

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(ENGLISH)

2. CHOOSE BETWEEN 6 MODES

- 1) Choose easily between 6 modes.
- 2) 12 min after the stopping exercise, your device will turn off automatically.

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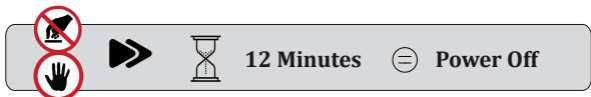
Recharge the battery.

Connect it to the USB cable.

Connect the USB cable to a computer or other USB port.

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(ENGLISH)

Assemble the host device on abdominal pad

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(ENGLISH)

Take off the adhesive part of the pad.

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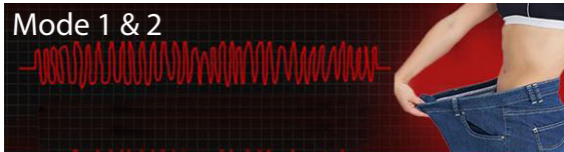
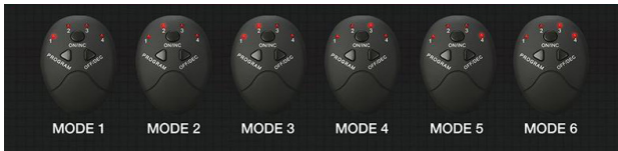
(ENGLISH)

Mode 1 & 2

In these two modes, muscles move intensely and burn fat quicker.

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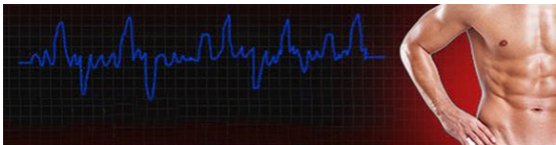
(ENGLISH)

Mode 3 & 4

Strong elastic movement and deep stimulation to help to promote muscle growth.

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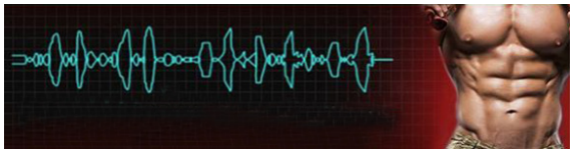
(ENGLISH)

Mode 5 & 6

These two modes are made to help you get a perfect figure.

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ENG (ENGLISH) This product must not be disposed of in municipal waste. To preserve the environment, this product must be recycled after its useful life as required by law.

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